

## Detergents & Cleaners

**Acidic  
Powder  
Liquid**

**Alkaline**  
Enzymatic Ready-to-use  
**Lens Cleaner**

Ultrasonics

**Neutral**  
Water Bath Conditioner  
**Liquid**  
Enzymatic Ready-to-use  
**Detergents**  
**Cleaners**

Detergents

**Concentrated Manual  
Cleaners**  
Concentrated  
**Ready-to-use**

Water Bath Conditioner

**Lens Cleaner**  
Concentrated

**Liquid**  
Acidic  
Enzymatic



460 Glennie Circle  
King of Prussia, PA 19406  
(800) 332-6647  
[www.deconlabs.com](http://www.deconlabs.com)

## MANUAL CLEANERS

### Contrad® 70

#### The Ultimate Labware Cleaner

Contrad® 70 is a liquid detergent concentrate that gets labware clean by soaking, scrubbing and ultrasonics. It rinses away totally, leaves zero residues. Phosphate-free, disposable straight to drain, Contrad® 70 is perfect for all critical cleaning, from labware to process equipment, ultrasonic cleaning and radioactive decontamination. The original and still the best... use it and see!



### Dri-Clean®

#### Powdered Detergent

Dri-Clean® is a powdered detergent for soak, manual, and ultrasonic cleaning applications. User friendly bottle provides easy dispensing and eliminates caking. Completely soluble, free-rinsing, and residue-free. Versatile formula removes a wide variety of tenacious contaminants.



### Neutrad®

#### Neutral pH Liquid Cleaner

Neutrad® is a liquid detergent concentrate for hand scrubbing your labware which offers high-performance cleaning with the safety of a near-neutral pH. Phosphate-free, totally rinsable, disposable straight to drain, Neutrad® is the perfect replacement for old-fashioned caustic soap powders.



### Citrad®

#### Acidic Liquid Detergent

Citrad® is a citric acid-based detergent concentrate ideal for removing salts, scale, and metal-oxide. Also used for metal brightening and pharmaceutical-process equipment. Completely phosphate-free, totally free-rinsing and biodegradable.



## AUTO WASHERS

### Dri-Contrad®

#### Powdered Autowasher Detergent

Dri-Contrad® is a powdered detergent that supercharges your automatic labware washer. Non-foaming, instantly soluble, rinsable, and chlorinated for biological control. Dri-Contrad is the only powder you need for your labware washer.



### Contrad® NF

#### Liquid Non-Foaming Autowasher Detergent

Contrad® NF is a non-foaming super-concentrated liquid detergent formulated for automatic labware washers. Phosphate-free, and totally free-rinsing, Contrad NF maximizes the efficiency of your labware machine.



### Contrex® AL

#### Alkaline Liquid Detergent

Contrex® AL is a concentrated alkaline liquid detergent formulated for use in manual, soak or ultrasonic cleaning of labware, scientific apparatus and process equipment. Contrex® AL is free-rinsing, phosphate-free and disposable straight to drain.



### Contrex® AP

#### Alkaline Powdered Detergent

Contrex® AP is a concentrated powdered labware detergent formulated for use in manual, soak or ultrasonic cleaning of labware, scientific apparatus and process equipment. Contrex® AP is completely soluble, free-rinsing and safe to dispose straight to drain. The user friendly bottle allows for easy dispensing and eliminates caking.



### DeSCAL®

#### Acidic Neutralizing Rinse (Liquid Descaler)

DeSCAL® is an inorganic acid-based detergent and descaler concentrate for removing hard water mineral deposits, rust and crusty build up from labware, water baths and stills. Low-foaming formula allows it to be used for automatic washing applications.



### Contrex® CA

#### Acidic Liquid Detergent

Contrex® CA is a concentrated acidic liquid detergent made for use in descaling, mineral removal, oxide removal or in applications where an acidic detergent is indicated. Formulated with organic acid, Contrex® CA can be used as a manual, soak or ultrasonic cleaner.



### Contrex® AW

#### Powdered Autowasher Detergent

Contrex® AW is a concentrated low-foaming powdered cleanser formulated for use in all automatic labware washing machines. Contrex AW is completely soluble, free-rinsing and biodegradable.



### Contrad® 100

#### Liquid CIP Cleaner

Contrad® 100 is a low-foaming detergent designed for clean-in-place systems and other high-pressure spray applications, including automatic washers. Alkaline formulation is phosphate-free, residue-free, and biodegradable.



See last page for additional information on Decon detergents.

## SPECIALTY CLEANERS

### ElimiNase®

#### RNase/DNase Remover

ElimiNase® Is a non-toxic, ready-to-use solution for instantly removing RNase, DNase, unwanted DNA, and radioactive contamination from labware and equipment. Safe and effective. Eliminates need for autoclaving or using DEPC.



### Contrex® EZ

#### Powdered Enzymatic Detergent

Contrex® EZ is a powdered detergent formulated for use as an enzymatic detergent for manual or soak cleaning of labware and equipment contaminated with proteins and fats. Contrex® EZ prevents corrosion, is free-rinsing and disposable straight to drain.



### Contrex® CF

#### Cation / Phosphate-Free Liquid Detergent

Contrex® CF is a concentrated liquid detergent formulated for use in cleaning of equipment and electronics that are sensitive to metal cations and phosphates. Contrex® CF is free-rinsing, phosphate-free and disposable straight to drain. Can be used for soak, scrub, sonication and auto washers.



### Optical Lens Cleaner

#### Optical Glass Surface Cleaner

Optical Lens Cleaner is excellent for removing oil and smudges from objectives or any optical surface. It gently cleans and will not affect a coated optic or soften the mounting that holds the objective in place. This lens cleaner is laboratory-safe and non-flammable. Regular cleaning helps extend the life of your objectives by removing surface debris that could cause scratching



### Jet Stream®

#### High Pressure Dust Remover

Jet Stream® High Pressure Dust Remover is moisture-free propelled air used to remove microscopic dust from negatives, microscope slides, lenses, balances and other laboratory equipment.



### BathBrite™


#### Water Conditioner for Water Baths

Bath Brite™ can be used to condition the water in water baths, incubators, and other temperature controlled lab equipment. 10 to 15 drops per gallon of water will keep water clear and fresh. Using Bath Brite™ will extend the time between water replacement.



## Detergents & Cleaners

additional information

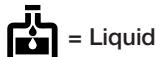

= Liquid



= Powder

DeSCAL<sup>®</sup>



Contrex<sup>®</sup> CA



Citrad<sup>®</sup>



Neutral<sup>®</sup>



Contrex<sup>®</sup> AL



Dri-Clean<sup>®</sup>



Contrex<sup>®</sup> AP



Contrad<sup>®</sup> 70



Dri-Contrad<sup>®</sup>



Contrad<sup>®</sup> NF



Contrex<sup>®</sup> AW



Contrad<sup>®</sup> 100



pH: ACIDIC

NEUTRAL

ALKALINE

SOAK, MANUAL SCRUB, ULTRASONIC CLEANING

ALKALINE

AUTO WASHERS

## PRODUCT GUIDE

PRODUCT	SIZE	DECON CAT. #	CASE SIZE
<b>Contrad 70® Ultimate Labware Liquid Detergent</b>	1 Liter 5 Liters 20 Liters 210 Liters	1002 1003 1004 1005	12 x 1 L 4 x 5 L — —
<b>Dri-Clean® Powdered Detergent</b>	2 kG 50 kG	2502 2501	6 x 2 kG —
<b>Neutrad® Neutral pH Liquid Cleaner</b>	950 mL 5 Liters 20 Liters	3003 3001 3002	12 x 950 mL 4 x 5 L —
<b>Citrad® Acidic Liquid Detergent</b>	3.8 Liters	4401	4 x 3.8 Liters
<b>Contrex® AL Alkaline Liquid Detergent</b>	1 qt. 3.8 Liters	5332 5301	12 x 1 qt. 4 x 3.8 Liters
<b>Contrex® AP Alkaline Powdered Detergent</b>	4 lbs. 25 lbs.	5204 5225	9 x 4 lbs. —
<b>DeSCAL® Acidic Neutralizing Rinse (Liquid Descaler)</b>	1 Liter 3.8 Liters 115 Liters	4001 4004 4030	12 x 1 L 4 x 3.8 Liters —
<b>Contrex® CA Acidic Liquid Detergent</b>	3.8 Liters	5701	4 x 3.8 Liters
<b>Dri-Contrad® Powdered Autowasher Detergent</b>	2 kG 5 kG 12.5 kG 50 kG	2002 2005 2007 2004	6 x 2 kG 4 x 5 kG — —
<b>Contrad® NF Liquid Non-foaming Autowasher Detergent</b>	3.8 Liters 10 Liters 115 Liters 210 Liters	6002 6001 6003 6055	4 x 3.8 Liters 2 x 10 L — —
<b>Contrex® AW Powdered Autowasher Detergent</b>	4 lbs. 25 lbs. 50 lbs.	5604 5625 5650	9 x 4 lbs. — —
<b>Contrad® 100 Liquid Clean-in-place Detergent</b>	3.8 Liters 20 Liters 115 Liters 210 Liters	1504 1503 1506 1505	4 x 3.8 Liters — — —
<b>ELIMINase® RNase/DNase Remover</b>	250 mL 950 mL	1101 1102	12 x 250 mL 12 x 950 mL
<b>Contrex® EZ Powdered Enzymatic Detergent</b>	4 lbs. 25 lbs.	5404 5425	9 x 4 lbs. —
<b>Contrex® CF Cation / Phosphate-Free Liquid Detergent</b>	3.8 Liters	5801	4 x 3.8 Liters
<b>Jet Stream® High Pressure Dust Remover</b>	11 oz.	3711	6 x 11 oz.
<b>Optical lens cleaner™ Optical Glass Surface Cleaner</b>	6 x 2 oz.	1100	8 x 6/pk
<b>Bath Brite™ Water Conditioner for Water baths</b>	8 oz.	3408	6 x 8 oz.