

# Microbiology & Histology

**Microbiology**  
**Histology**

**Cytology Fixative**  
**Clearing Agents** **Odor Control**

**Loops**  
Needles

**Nichrome**

**Formalin Neutralizer**

**Freezing Media**

**Loops**

**Needles**

**Streakers**

**Nichrome**  
Odor Control

**Calibrated**

**Formalin Neutralizer**

**Calibrated**  
**Non-calibrated**  
Clearing Agents

**Streakers**

**Freezing Media**  
**Needles**  
Cytology Fixative



460 Glennie Circle  
King of Prussia, PA 19406  
(800) 332-6647  
[www.deconlabs.com](http://www.deconlabs.com)



**CitriSolv™**  
*Solvent and Clearing Agent*

CitriSolv™ is a d-Limonene-based solvent and clearing agent that can be used as a safe alternative to xylene and ethyl acetate. It is soluble with absolute alcohol and will dissolve paraffin, glues, adhesives, and mounting media. CitriSolv provides for better tissue sections – it will not harden or shrink tissue. It dries slowly with no residue and enhances some stains. Use in manual or automatic tissue processor staining procedures in histology/cytology applications that require xylene.

Item	Standard Size	Alternate Size	Decon Catalog #
CitriSolv™	1 G (3.8L)	4 x 1 G (3.8L)	1601



**CitriSolv™ Hybrid**  
*Solvent and Clearing Agent*

CitriSolv™ Hybrid is a mixture of d-limonene and hydrocarbon material allowing for the reduction of orange aroma and an increase in flash point. Not classified as flammable.

Item	Standard Size	Alternate Size	Decon Catalog #
<del>CitriSolv™ Hybrid</del>	<del>1 G (3.8L)</del>	<del>4 x 1 G (3.8L)</del>	<del>1601H</del>



**HydeAway®**

HydeAway® aldehyde decontamination management products effectively convert used aldehydes into non-hazardous material to ensure safety and to facilitate disposal. HydeAway® offers the highest neutralization ratio available and can be used for formaldehydes, formalin and glutaraldehydes. Crystalline Hyde-Away® turns aldehydes into a solid polymer at a 1:30 ratio and the liquid version turns aldehyde into a fluid suitable for sewer system disposal at a 1:15 ratio. Tested via EPA methods and proven safe for drain or trash disposal. Treating used aldehyde with Hydeaway™, prior to drain disposal, is better for the environment and makes for a safer workplace.


Item	Standard Size	Alternate Size	Decon Catalog #
Hyde-Away® (Crystals)	30 oz.	12 x 30 oz.	2224
Hyde-Away® (Liquid)	1 G (3.8L)	4 x 1 G	2201



**Flash Freeze™**

Flash Freeze™ is moisture-free propelled cold gas used to rapidly cool or freeze samples in clinical or scientific laboratories.

Item	Standard Size	Alternate Size	Decon Catalog #
Flash Freeze™	11 oz.	6 x 11 oz.	3511

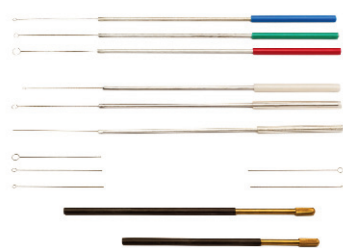

### Inoculating Loops, Streakers and Needles

Made of high quality, twisted nickel chromium wire, LeLoops and LeStreaks are re-usable—making them more economical than disposable plastic loops. Loops slide smoothly across the agar surface

without tearing or gouging. The nichrome wire cools rapidly after flaming, saving time. Loops are available calibrated or non-calibrated. All calibrated loops are subjected to Evan's Blue dye test.

LeLoops come without holders, which are available separately. Holders are made of nickel-plated brass with well insulated handles.

LeStreak and LeStab are permanently attached to lightweight, insulated aluminum holders perfectly balanced for easy manipulation.

<b>LeLoops®</b>				
Item	Loop ID	Std. Size	Alt. Size	Decon Cat #
Calibrated Loop 10 µL		5/pk	5pk/cs	MP190-25
Calibrated Loop 1 µL		5/pk	5pk/cs	MP195-25
Small Loop	2.26mm	25/pk	—	MP197-25
Medium Loop	2.91mm	25/pk	—	MP199-25
Large Loop	3.57mm	25/pk	—	MP198-25
<b>LeStreak®</b>				
Small Streaker Blue	2.26mm	5/pk	—	MP184-5
Medium Streaker Green	2.91mm	5/pk	—	MP182-5
Large Streaker Red	3.57mm	5/pk	—	MP180-5
Calibrated Streaker 1 µL		5/pk	5 x 5/pk	MP185-5
Calibrated Streaker 10 µL		5/pk	5 x 5/pk	MP187-5
<b>LeStab® Needles</b>				
LeStab Micro Needle		5/pk	—	MP186-5
<b>LeLoop®/ Loop Holders</b>				
Brass Loop Holder 6"		1 each	10/pk	MP6-10
Brass Loop Holder 8"		1 each	10/pk	MP8-10



### Odorgone® & The Pill™

*Autoclave Deodorizer*

Heat and moisture from the autoclave release the contents of the capsule to overcome unpleasant smells generated from heating microbiological waste. Simply place

one capsule in each load to absorb offensive odors leaving the laboratory smelling fresh and clean. Odorgone comes in two pleasant scents; Lemon and Apple. The PILL is fresh scented.

Item	Standard Size	Alternate Size	Decon catalog #
OdorGone—Apple	25 Capsules	8 x 25 Capsules	APP25
OdorGone—Apple	100 Capsules	6 x 100 Capsules	APP100
OdorGone—Lemon	25 Capsules	8 x 25 Capsules	LEM25
OdorGone—Lemon	100 Capsules	6 x 100 Capsules	LEM100
The PILL	100 Capsules	20 x 100 Capsules	5100
The PILL	500 Capsules	6 x 500 Capsules	5500



### SlideRite™ Cytology Fixative

SlideRite™ is a highly concentrated, water-soluble spray that forms a lasting protective coating over all types of cytological smears. Many cytological specimens are fragile and degenerate quickly, therefore, an immediate fixation of the material to be studied is essential.

Item	Standard Size	Alternate Size	Decon Catalog #
SlideRite	1.5 oz.	12 x 1.5 oz.	1801